

### LED Backlight

Led backlight for viscometer baths TV2000 - TV4000 - TV7000 - TV2500



Backpanel

Stand alone



- Low energy consumption (6 Watt)**
- Long lifetime (LED's)**
- Wide range input**
- Easy to mount**
- Excellent contrast**

Item	Unit	All models
Operating voltage	[V]	85 - 264
Operating frequency	[Hz]	47 - 63
Weight	[kg]	2,5
Ambient condition	[°C]	18 .. 28
CE	All models conform to CE regulation	

Z41 stand alone		00T0909
Weight	[kg]	2,5
Power	[Watt]	6
Dim LED. LxW	[mm]	244 x244

Z41 backpanel		00T0908
Weight	[kg]	2
Power	[Watt]	6
Dim LED. LxW	[mm]	244 x244

Z71 backpanel		00T0907
Weight	[kg]	3,5
Power	[Watt]	12
Dim LED. LxW	[mm]	2 x 244 x244

#### General

To improve visibility inside the bath Tamson offers a LED based backlight. There is a choice between a lamp with adjustable height (stand alone) or a specific backlight panel which can be mounted directly onto the bath.

#### LED technology




The used LED technology offers several advantages. The system has a long life expectancy of over 20.000hrs. Power consumption is limited drastically to just several watts. Optimal light distribution minimises dark or light spots

#### Installation

All models offer a connection via a standard mains (European) power plug connection and can be connected to the mains voltage. The backpanel can be powered from the additional powerplug located on the back of our TV models. This plug standard is included in the TV.


### LED Backlight

Led backlight for viscometer baths TV2000 - TV4000 - TV7000 - TV2500

Item	Picture	Article number and description
00T0909		Stand Alone LED Panel Wide range input (85 .. 230V AC, 50-60Hz) For use with TV2000, TV4000, TV2500 or stand alone
00T0908		LED Panel "Back mount" wide range input For use with TV2000 and TV4000 only
00T0907		LED Panel "Back mount" wide range input For use with TV7000 only

Mounting the backpanel:

- 1: LED panel
- 2: Power plug on panel (euro) female
- 3: Power plug on bath (euro) female
- 4: Powercord, (euro) male/female

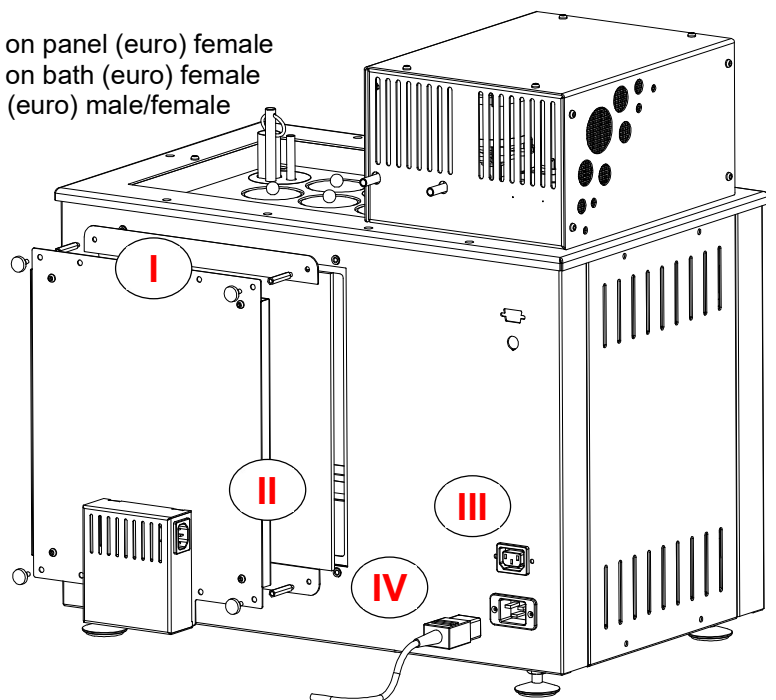

Z41-71